

HARMAR SL600 PINNACLE STAIR LIFT



The World's Most Compact Stair Lift

The Pinnacle is the slimmest and most light-weight stair lift on the market. This compact stair lift folds to take just 10.5" of staircase width, so everyone in the house can breeze by it with ease. Enjoy the smoothest, safest, and most energy-efficient ride by adding this low-profile stair lift to your home.






The patented worm drive requires no

greases or lubricants.




SAFE Safety sensors on chassis and footrest stop the lift instantly at obstacles.

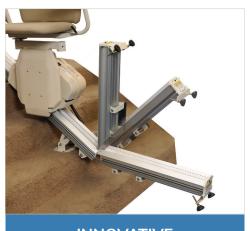

INNOVATIVE Industry exclusive, optional fully automatic fold rail design eliminates obstructions at the bottom of the stairway.



HARMAR SL600 PINNACLE STAIR LIFT

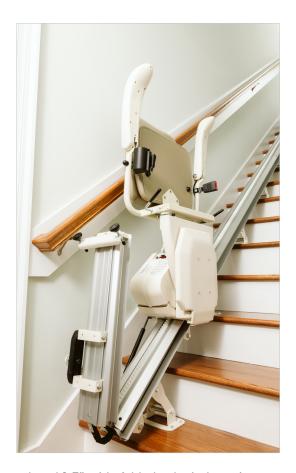

- Just 10.5" wide folded the industry's narrowest
- Patented helical worm drive is smooth & energy-efficient
- DC drive even runs during power interruptions
- Clean installation and maintenance
- Powerful enough to travel up to a 52° incline
- Patented Folding Rail option eliminates obstruction at stair base (optional)
- Includes 2 wireless remotes
- ETL Listed 3148125, Commercially rated

SPECIFICATIONS	SL600
Maximum Weight Capacity	350 lbs.
Overspeed Governor	Standard, Commercially Rated
Track (rail) Type	Clear Anodized Aluminum
Travel	15'6" Standard; 70' Maximum
Average Number of Return Trips per Charge	Up to 60
Control in Armrest	Yes (left or right hand operation)
Lift Mounts to Steps or Wall	Steps
Minimum Folded Width	10.5"; 13.6" with Folding Rail
Minimum Footrest Height	3.5"
Clear Distance Between Armrests	19"
Floor to Seat Height Low / High	20" / 24.5"
Minimum Wall to Stair Side of Rail	6"
Seat Depth	16"
Backrest Height From Top of Seat	18"
Electrical Requirements	120-240 V AC
Operation Power	24V DC Battery
Speed	20 fpm
Incline Range	27° - 52°; Inclines above 45° may require local jurisdiction approval
Drive System	Helical Worm Gear with Nylon Polymer Rack
Folding Rail Option	Automatic and Mechanical
Safety Design	Complies with ASME A18.1, CAN/CSA B44.1, ASME A17.5
Safety Devices	 Centrifugal Overspeed Brake Bi-directional Footrest Obstruction Sensor Swivel Seat Cut-off Switch Final Limits/Seat Belt Chassis Obstruction Sensors Optional Key Switch
Warranty	Limited Lifetime Warranty 1-Year: Batteries



1 Year Limited Warranty on New Equipment Installs

An independently owned and operated franchise business.



⊘ Vertical Platform Lifts **⊘** Pool Lifts

⊘ Commercial

Patient Lifts



Auto Lifts

Bath Modifications

And More