



The Alzaro is the perfect choice when a strong, sturdy floor lift is required. The Alzaro features a curved boom which allows for a greater lifting height without the user getting too close to the lift. With an optional balancing cradle or a loop spreader bar, the Alzaro ensures all lifting requirements can be met.

Unique Features

The Alzaro features a single frame design that ensures the lift is extremely strong, highly stable and very reliable. The legs open and close electronically and the battery can be plugged in for external charging. It is ideal for hands-free lifting and transferring of patients up to 160kg/352.7lb. With a hexagon-shaped push handle, and an overall width of 68.7cm/27", the Alzaro can easily be moved by one caregiver for use in multiple rooms. The high quality casters ensure smooth rolling and the rear casters lock to prevent accidental movement.

Optional Accessories

The Alzaro comes standard with a loop spreader bar but can easily be fitted with a

self-balancing cradle for clip slings (manual or powered), or a weight scale.

Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B
117 43 Stockholm
Phone: +46 8 665 35 00
Fax: +46 8 665 35 10
info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street
Port Melbourne, VIC 3207
Phone: +61 1300 132 736
Fax: +61 (03) 9645 4954
info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road
Ottawa, ON K2E 7K1
Phone: +613.723.6734
Fax: +613.723.1058
info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124
8076 PA Vierhouten
Phone: +31 577 412 171
Fax: +31 577 412 170
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K
Austin, TX 78754
Phone: +512.476.7199
Fax: +512.476.7190
info.us@humancaregroup.com

Accessories




Spreader Bar

Standard with lift. Fits most loop slings. Upholstered cover may be used (Part no. 50866)



Self-Balancing Cradle

Manual or electric. Fits clip slings. Must be ordered with lift at time of purchase.

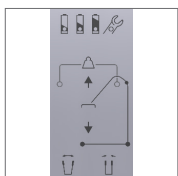

Hand Control

Four functions. Included with lift. Lift arm up/down, electrical leg spreading.



Wheels

Rear wheels feature kick-brake. Diameter 125mm.

Lifting Situations




Use to/from a bed

Use to raise/lower a patient from/to a bed.



Use to/from a chair

Use to raise/lower a patient from/to a chair.

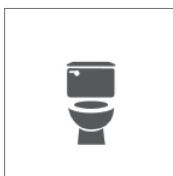

Use to/from the toilet


Use to raise/lower a patient from/to a toilet.

Features




Curved Boom

Allows for a greater reaching height without getting too close to the user.

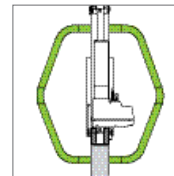

Hexagonal Push Handle

Comfortable width makes it easy to move lift from room to room.

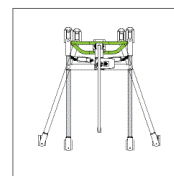

Spreadable Base

Electrical adjustment from handset or control box.



Tested to all relevant standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Alzaro			
Part No.	11000490	IP Rating	IPX4
Width	68.7cm / 27"	Turning Radius	148cm / 58.3"
Length	118.3cm / 46.6"	Front Wheels	10cm/4"
		Rear Wheels with Brake	12.5cm/5"
Height	Arm Lowered: 86.3cm / 34" Arm Raised ¹ : 195.3cm / 76.9"	Battery	Pb battery 24V DC/2.9Ah, 2x 12V/2.9Ah Li-ion battery
Height from Floor	11.3cm / 4.4"	Weight	45kg / 99.2lb
Adjustable Leg (Min - Max)	Outside: 67.9cm - 119.6cm / 26.7" - 47.1"	Weight Capacity	160kg / 352.7lb
Sound Level	48dB (loaded) while lifting 52dB (loaded) while spreading legs		

¹Measurement based on arm mast at highest setting